



Modular membrane adsorber based system for protein purification

Öznur Kökpınar¹, Fabienne Anton¹, Daniela Harkensee¹, Cornelia Kasper¹, Robert Zeidler², Oscar-Werner Reif³, Roland Ulber⁴ and Thomas Scheper¹

¹Institut für Technische Chemie, Callinstr. 3, 30167 Hannover, Germany

²Vivascience AG, Feodor-Lynen-Str. 21, 30625 Hannover, Germany

³Sartorius AG, Weender Landstrasse 94-108, 37075 Göttingen, Germany

⁴AG Bioverfahrenstechnik, Fachbereich Maschinenbau und Verfahrenstechnik, Technische Universität Kaiserslautern, Germany

Abstract

The aim of the project was the development of a modular construction system concept based on membrane adsorber technology for protein purification. The approach enables a high throughput screening for the determination of the optimal downstream conditions of several bioproducts simultaneously. An easy scale up of the purification results from microtiter to industrial scale can be performed.

HTDSS

In order to screen different target proteins and to devise the most efficient strategy for their fast and direct separation and purification, we developed a modular construction system concept based on membrane adsorber technology. Figure 1 exemplifies this concept of the so called high throughput downstream screening system (HTDSS).

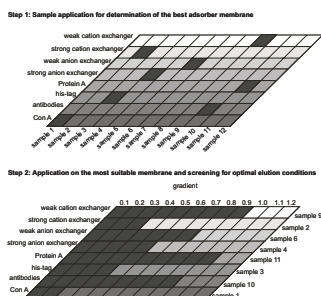


Figure 1. Modular concept for a screening experiment

Easy Scale-up

For the implementation of the HTDSS 8-strip devices were designed. Each 8-strip represents a certain membrane type (e.g. cation exchanger, metal chelate) and can be arranged in the 96-well format and processed by centrifugation.

After the determination of the optimal purification strategies for the target proteins scale-up experiments can be performed using membrane modules in various formats (Figure 2).

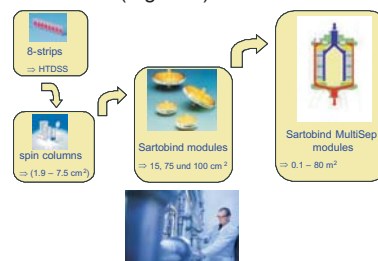


Figure 2. Membrane adsorber in various formats

Separation of model proteins


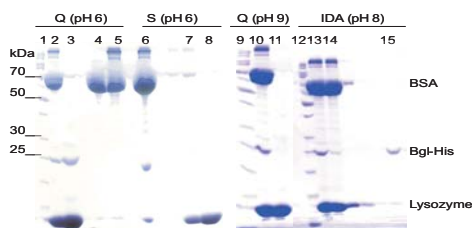


Figure 3. Separation of BSA, lysozyme and Bgl-His on Q, S and IDA membrane adsorber.

The efficiency of the used adsorber membranes were tested by the separation of a BSA/lysozyme/Bgl-His. The different pI values and the distinct affinities for the metal chelate allow the separation of BSA, lysozyme and Bgl-His on ion exchange and IDA membrane adsorbers. The isolation of BSA was carried out at pH 6 on Q membrane (Lane 4). S membrane retained lysozyme at this pH value (Lane 7). Lysozyme was also fractionated at pH 9 on Q membrane, where BSA and Bgl were completely adsorbed on the membrane and pure lysozyme passed through (Lane 11). The isolation of Bgl-His from the mixture succeeded only with IDA membrane. Bgl-His was purified on IDA membrane loaded with Cu^{2+} -ions (Lane 15).

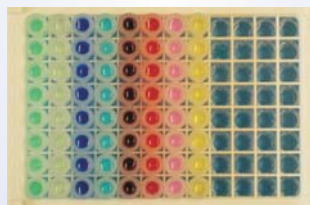

Purification of IgG₁

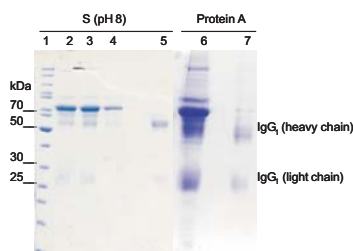


Figure 4. Purification of IgG₁ from cell culture supernatant using S and Protein A membrane.

As an other example for the successful application of HTDSS the purification of IgG₁ from cell culture supernatant is presented. The antibody could be purified on S membrane at pH 8 (Lane 5) and by applying on Protein A membrane (Lane 7).

Purification of GFP-His

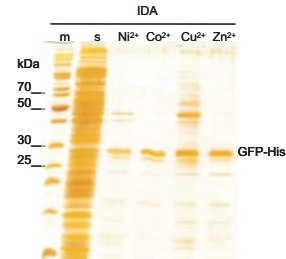


Figure 5. Purification of GFP-His using metal chelate membrane.

The purification of recombinant proteins engineered with polyhistidine, (His)₆ tags from crude fermentation broths and cell lysates using Immobilized Metal-Affinity Chromatography (IMAC) provides a one-step isolation of proteins at over 90 % purity. 8-strips allow a screening of different metal ions and matrices simultaneously. To obtain pure GFP-His, the 8-strips were loaded with Ni^{2+} , Co^{2+} , Cu^{2+} and Zn^{2+} solutions and the His-tagged protein was separated from cell lysate by applying on the charged membranes. According to the analysis by SDS-PAGE, Co^{2+} -ions are most suitable for the isolation of GFP-His.

Conclusion

The applicability of the developed membrane adsorber based devices in HTDSS, the selectivity of the adsorber membranes and the effect of loading-elution conditions were investigated and successfully applied for the purification of different proteins from model protein mixture, fermentation broth and cell culture supernatant.

In summary the HTDSS facilitates a fast and parallel determination of the optimal purification conditions of the target proteins.

Acknowledgment

This research project was supported by Deutsche Bundesstiftung Umwelt (DBU, AZ 13099). We thank PD Karl Friehs, AG Fermentationstechnik, University of Bielefeld for providing of the bacterial cultures.